


SofCare® DuoGard® 2-Layer Reactive Overlay

Convenient. Complete. Effective.

Stryker's SofCare family of products allows for a convenient, complete and effective way to assist pressure redistribution while improving hospital outcomes by helping to prevent facility acquired pressure ulcers.

With its 2-layer design, the DuoGard Overlay provides economical protection by lowering tissue interface pressure below capillary closure to help prevent and treat pressure ulcers. Its soft, durable vinyl construction easily wipes clean and is simple to use. Its convenient hand wells facilitate patient transfer and repositioning.

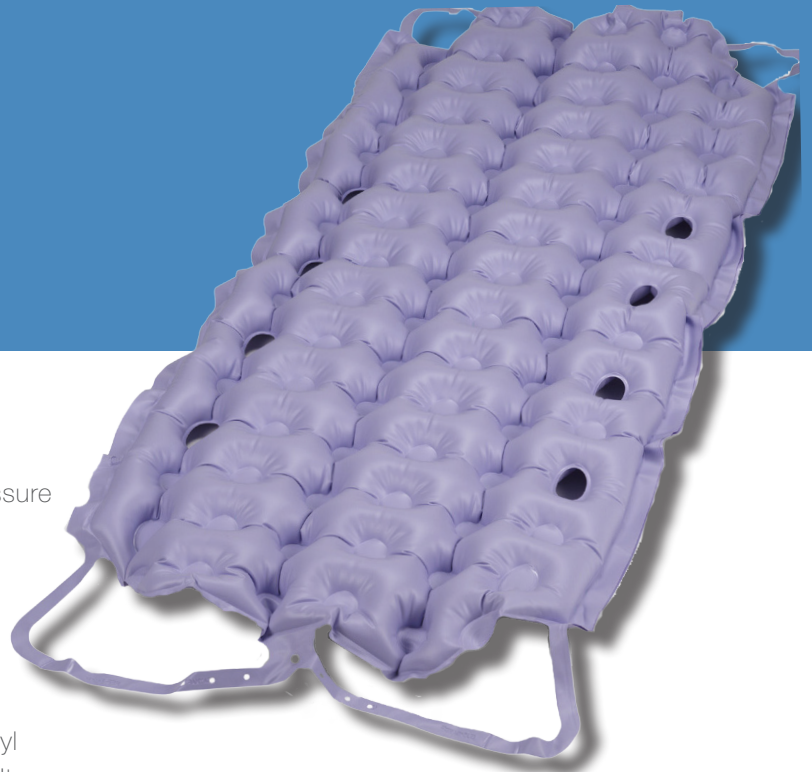

Inflation Options




The **Companion Pump 2.0** is a constant low pressure powered air pump intended to provide continuous air pressure to the overlay to assist in the treatment and prevention of pressure ulcers.



The **SofCare Manual Hand Pump** is a non-powered inflation device that can go home with the patient.


Air-Channelling Technology

SofCare air-channelling technology redistributes tissue interface pressure to aid in prevention and treatment of pressure ulcers. The channel design allows for air circulation helping to manage patient microclimate.

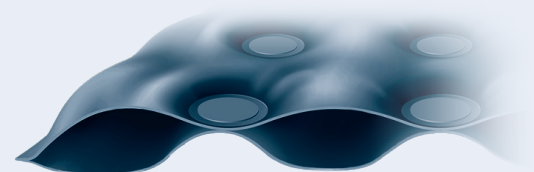

SofCare Air-Channelling Technology

SofCare[®] DuoGard[®] Overlay

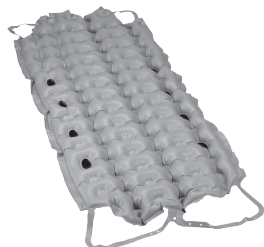

SofCare DuoGard Features

- Channeled tubular support
- 2-layer design
- Vinyl material
- Hand wells/Transfer handles
- Corner straps
- 2 inflation options

Specifications

	DuoGard Overlay
Model Number	2780-300-000
Inflated Length	76" (193 cm)
Inflated Width	33" (83.8 cm)
Weight	4.5 lbs (2.04 kg)
Weight Limit	400 lbs (181.4 kg)
Packaging	6/carton
Materials	Vinyl
Warranty	90 days

Pump Options

2781100000	Companion Pump 2.0
HCP700	SofCare Manual Hand Pump

Single Patient Use Only

Ask about our other SofCare 3-layer air-channeling technology products. Stryker reserves the right to change specifications without notice.

3800 E. Centre Ave.
Portage, MI 49002 U.S.A.
t: 269 389 2100 f: 269 329 2311
toll free: 800 787 9537

www.stryker.com