

PositionPRO™

Patient Repositioning

Non-integrated, low air-loss, pressure redistribution surface helps prevent and treat skin ulcers

Infinite positioning options up to 40 degrees with a hold time of up to 120 minutes

Turn-assist and maximum inflate capabilities for patient positioning controlled by intuitive user pendant


Surface rotates while bed is flat or with the head end raised up to 45 degrees



Attached pendant conveniently clips anywhere on the bed. Intuitive control provides one-touch patient positioning.

Infinite Patient Positioning Options

PositionPRO is the only product in its class that can provide full position assistance while providing powered pressure redistribution, with an industry-leading 40 degree turn and 120-minute hold. This requires much less force to facilitate a turn, resulting in reduced back injuries and increased efficiency for nurses, and the reduction of pressure ulcers for patients.

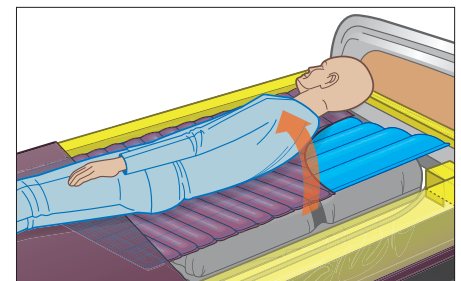

Industry-leading 40-degree turn with 120-minute hold time.

PositionPRO™

Patient Repositioning

stryker®



Medical

3800 E. Centre Ave.
Portage, MI 49002
t: 269 329 2100 f: 269 329 2311
toll free: 800 787 9537

www.stryker.com

Standard Features

- Turn-assist capability with infinite positioning
- Low air-loss
- Maximum inflate
- Three layers of pressure-redistributing surface
- Dartex® cover
- Fluid and bacterial resistant cover
- 16CFR1633 and 1632 Compliant
- CAL TB 603 Compliant
- Firm foam perimeter

Specifications

• Model Number	2920
• Overall Length	84" (213 cm)
• Overall Width	35" (89 cm)
• Thickness	6.5" (16.51 cm)
• Weight Capacity	500 lbs (228 kg)

*Stryker reserves the right to change specifications without notice.
Support surfaces are manufactured for Stryker Medical.
Dartex® is a registered trademark of Dartex Coatings, Ltd.*

Warranty

- Two-year on all components; includes parts, labor and travel

Electronics

• Agency Approvals	US UL
• Volts	120 VAC
• Ampere Ratings	2.5 A
• Current Leakage	< 300 microamperes
• Frequency	60 Hz
• Hospital Grade Plug	Yes
• EMC Compatible	Yes

Optional Features

- Nylon cover