

Sof.Care[®] DuoGard[®] Overlay



Convenient. Complete. Effective.

2-layer reactive overlay

Our Sof.Care family of products allows for a convenient, complete and effective way to assist pressure redistribution while improving hospital outcomes by helping to prevent facility acquired pressure injuries.

With its 2-layer design, the DuoGard Overlay provides economical protection by lowering tissue interface pressure below capillary closure to help prevent and treat pressure injuries. Its soft, durable vinyl construction easily wipes clean and is simple to use. Its convenient hand wells facilitate patient transfer and repositioning.




The Companion Pump 2.0 is a constant low pressure powered air pump intended to provide continuous air pressure to the overlay to assist in the treatment and prevention of pressure injuries.



The Sof.Care Manual Hand Pump is a non-powered inflation device that can go home with the patient.


Air-channelling technology

Sof.Care air-channelling technology redistributes tissue interface pressure to aid in prevention and treatment of pressure injuries. The channel design allows for air circulation helping to manage patient microclimate.

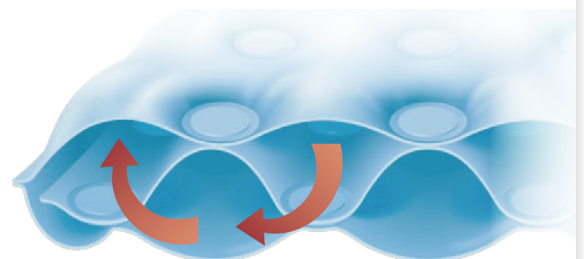

Sof.Care air-channelling technology

Sof.Care DuoGard Overlay

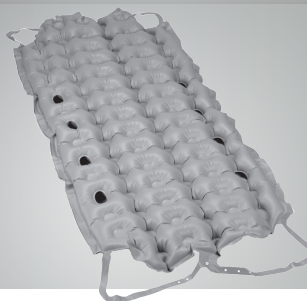

Sof.Care DuoGard features

- Channeled tubular support
- 2-layer design
- Vinyl material
- Hand wells/transfer handles
- Corner straps
- 2 inflation options

Specifications	DuoGard Overlay
Model number	2780-300-000
Inflated length	76" (193 cm)
Inflated width	33" (83.8 cm)
Weight	4.5 lb (2 kg)
Safe working load*	400 lb (181.4 kg)
Packaged	6/carton
Materials	Vinyl
Warranty	90 days

Pump options	Description
2781100000	Companion pump 2.0
HCP700	Sof.Care manual hand pump

*Therapeutic working load is 100 lb. (45.3 kg) to 350 lb. (158.8 kg).

Single patient use only.

Ask about our other SofCare 3-layer air-channeling technology products.

Stryker reserves the right to change specifications without notice.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: DuoGard, Sof.Care, Stryker. All other trademarks are trademarks of their respective owners or holder.

3800 E. Centre Avenue
Portage, MI 49002 USA
t: 269 329 2100
toll free: 800 327 0770

patientcare.stryker.com